


ZANDAIR™ AIR PURIFIERS

AIR PURIFICATION FOR IVF LABS

PROTECT YOUR LAB WITH EFFECTIVE AIR FILTRATION

OVERVIEW

No matter how well you maintain your laboratory, your work area can still be filled with air pollutants, volatile organic compounds (VOCs), aldehydes, carbon monoxide, viruses, and bacteria. Through use of the internal activated charcoal and high efficiency filters and Photocatalytic Oxidation Chamber, the ZANDAIR Air Purifiers help to reduce microparticles in the air and neutralize VOCs.

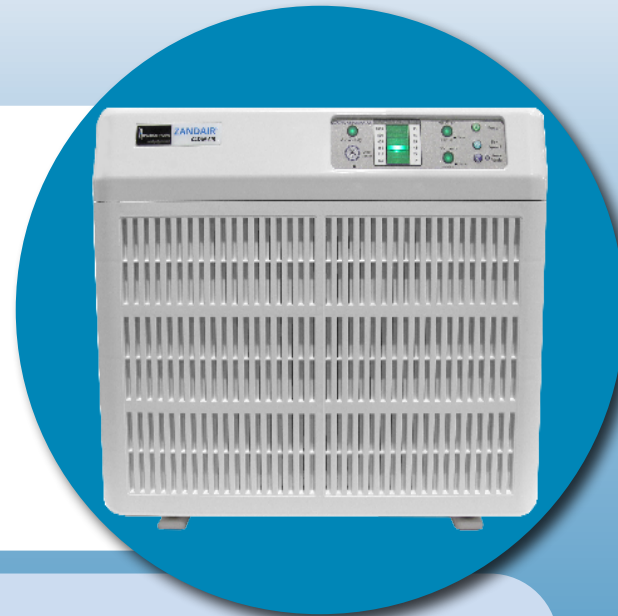

ZANDAIR 100C Portable Air Purifier


The ZANDAIR 100C unit is encased in a high impact ABS plastic with no off-gassing and is light weight and portable. Electronic sensors monitor air quality and automatically increase the performance of the air purification system. Warning lights alert staff to the presence of toxic chemicals and fumes well before they reach dangerous levels

ZANDAIR PCOC3 In-Duct Air Purification System




An encased unit to be fitted into the HVAC/AHU unit or the air ducting to provide an airtight flow of the airstream through the unit and has the option for a removable MERV-13 and Activated Charcoal (ACF) Filter. This model consists of three photocatalytic grids to provide additional dwelling time for maximum breakup of the molecular structure of the VOCs passing in the airstream.

ZANDAIR Features

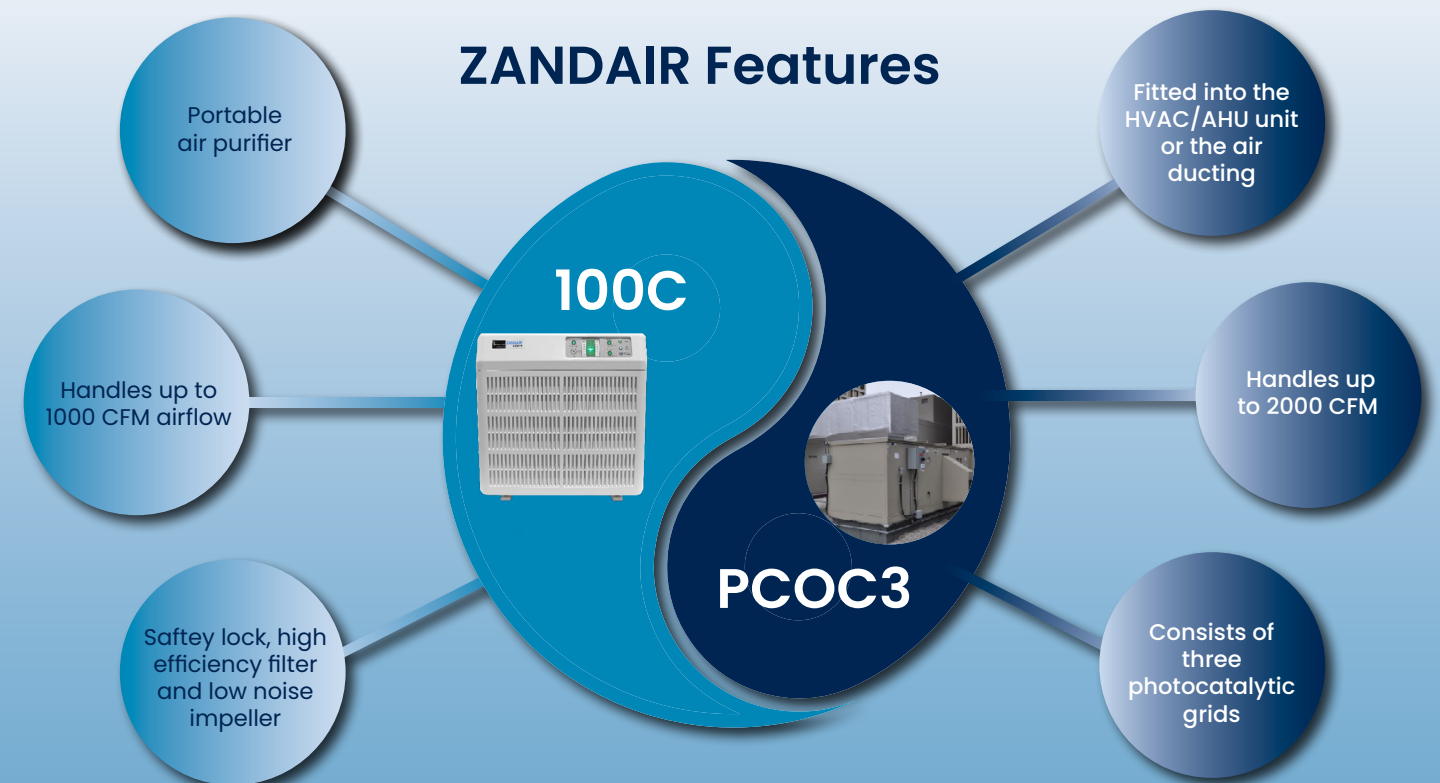

Photocatalytic Oxidation Process Overview

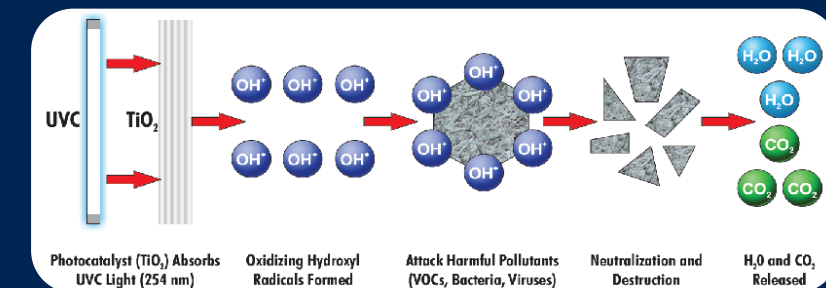

100C

Technical Specifications

PCOC3

- Power: 120V-60Hz/220V-50 Hz
- Max Air Flow: 265 CFM/7.5 m³ per minute
Application: 1000 ft³/29 m³
- Dimensions (w x h x d): 21.5" x 18" x 8" / 55cm x 46cm x 20cm, Weight: 23 lbs/10.43kg
- UVC Range: 254 nm
Service: Replace UVC Lamps-1 year/ Filters-6 months

- Power: 110V/60Hz / 220V/50Hz
- Max Air Flow: 2000 CFM/ 56.6 m³/minute
Application: 2000 ft³/57 m³
- Dimensions outside (w x l x h): 658 mm x 785mm x 585mm/25 7/8" x 30 7/8" x 23 7/8"
Weight: 92 lbs/42kg
- UVC Range: 253.7 nm

Contact our team to transform your lab today!

Hamilton Thorne Inc.
100 Cummings Center, Suite 465E
Beverly, MA 01915
USA
978-921-2050
www.hamiltonthorne.com
info@hamiltonthorne.com



Visit our website!